


Cheat Sheet

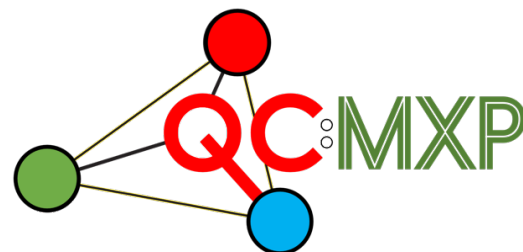

Prefilter

Option	Action
% allowable missing QCs <=	Metabolite feature is kept if there are less than the specified % missing QC samples.
↓ Every Batch	% missing QC calculated per batch and peak kept if every batch has <= specified % missing QCs
↓ Any Batch	% missing QC calculated per batch and peak kept if any batch has <= specified % missing QCs
↓ Complete	% missing QC calculated across all batches
<input checked="" type="checkbox"/> < blank = missing	consider any QC sample with signal < the blank threshold to be missing

Configuration file

Name	Options	Explanation
LogTransform	"true"/"false"	Perform Log10 transformation on data before correction (to ensure normal distribution)
RemoveZeros	"true"/"false"	There is no such thing as zero. Samples is either detected (value) or not detected (missing). Zeros mess up statistics.
OutlierMethod	"None", "Percentile", "Linear", "Quadratic", "Cubic"	The 'linear'/quadratic/cubic' options fit a simple polynomial function to the QC data using robust least squares regression. Points outside the population confidence interval are deemed outliers. The 'percentile' option implements the standard non-parametric < Q1 - 1.5 IQR or > Q3 + 1.5 IQR outlier detection method on the QC samples.
OutlierCI	value between 0.9 and 0.99	Set the Confidence interval for the polynomial methods above. 90%-99%
OutlierPostHoc	"Ignore", "MPV", "NaN"	For the corrected data. Either ignore outliers or replace with the QC median peak value (MPV) or remove (replace with "not a number" NaN).
IntraBatchMode	"Mean", "Linear", "Spline", "Sample"	Three correction modes. "Spline" is the default QCRSC algorithm that requires optimisation of the smoothing parameter. "Linear" is a simple Robust (bisquare) linear regression based on the QC values & requires no smoothing optimisation. "Mean" equalises the QC mean across batches & ignores within batch systematic change. "Sample" ignores the QCs and corrects using linear regression based on the samples labelled 'Sample' ('Sample' option is not recommended).
InterBatchMode	"QC", "Reference", "Sample"	'QC' = both the intra- and inter-batch based on the QC samples; 'Reference' = intra-batch correction based on the QCs but the inter-batch correction uses the samples labelled as 'Reference'. 'Sample' = both the intra-batch & inter-batch correction uses the samples labelled as 'Sample' ('Sample' option is not recommended).
QCRSCgammaRange	"x:y:z"	The range of values (x to z) to search to determine to optimal correction curve (in increments of y). A value >=4 results in a curve equivalent to a linear regression. A value < 0 results in a highly nonlinear curve.
QCRSCcvMethod	"3-Fold", "5-Fold", "7-Fold", "Leaveout"	Type of cross-validation used for optimising the smoothing parameter value.
QCRSCmcReps	integer	The number of Monte Carlo (random) resamples of the k-fold cross validation. The resulting cvMSE the mean of the generated set of cvMSEs.
CorrectionType	"Subtract", "Divide"	Subtract or divide the correction curve from the raw data. Optional depending on whether you believe the bias to be additive or multiplicative. Recommendation: If multiplicative & Log10 transformation is used then Subtract.
BlankRatioMethod	"QC", "Median", "Percentile"	BlankRatio = 100*max(BlankValue)/SampleReferenceValue. 'QC' sets the SampleReferenceValue as the median QC value; 'Median' sets the SampleReferenceValue as the median Sample value; 'Percentile' sets BlankRatio = % of Samples < max(BlankValue)*RelativeLOD (RelativeLOD is a user defined constant - default 1.5). In this instance missing values are considered Blank values
RelativeLOD	value	The RelativeLOD is a multiplier (relative to the Blank) used to set a value below which a peak is considered background noise. i.e., "estimated LOD = max(Blank) x relativeLOD"
StatsParametric	"true"/"false"	Choose whether summary statistics for each corrected feature are calculated using parametric (mean, standard deviation etc) or non-parametric methods (median, median absolute deviation etc).
ParallelProcess	"true"/"false"	Switch on the parallel processor to use multiple processor cores to speed up the QCRSC engine.

Clean & Explore

Name	Explanation
Missing Filter:	Remove peaks where there are low number of total samples.
Every-Batch/Any-Batch/Complete	Should the filter bank calculate its peak-wise statistics across all batches (complete) or calculate for each batch individually and if 'Every Batch' then every batch must pass the % missing threshold for a given peak to be retained, else if 'Any Batch' then only one batch as to pass the threshold for the peak to be kept.
% missing thesh =	% allowed missing Samples. e.g. if 20% then a given peak is removed if number of missing samples labeled 'Sample' > 20%.
Before/After	Compare the performance statistics/visualisations before or after correction. Use this to convince yourself of the utility of the correction algorithm.
Filter Bank:	
ON/OFF	Apply the dial setting to the peak selection filter.
Mode	leave on 'Complete' unless performing large scale studies over many batches. For a given peak: 'Max qcRSD' deferes to the using the stats from the worst batch (batch with highest qcRSD) for all filter settings. 'minQC' deferes to the best batch (batch with lowest qcRSD) for all filter settings. 'Median qcRSD' deferes to the most representative batch for all filter settings. 'Complete' calculates the stats across all batches.
PCA preprocessing:	
log10 Transform	apply a log10 transformation to all data.
Autoscale / Pareto scale	Preform autoscaling (mean centre then divide by standard deviation) or Pareto Scaling (mean centre then divide by sqrt(standard deviation) to each individual peak. Scaling is always performed after the transformation.
PCA Missing value imputation:	
KNN column	KNN missing value imputation replacing with the nearest metabolite feature
KNN row	KNN missing value imputation replacing with the nearest sample
k =	replaces missing values with a weighted mean of the k nearest-neighbor columns/rows. KNN imputation always performed after transformation & scaling.
blank/20%min	replaces missing values with the maximum blank value for that peak or if no blanks detected 20% of the lowest value. blank/20%min is always performed before transformation or scaling.
PCA projection:	
Project in QC samples	PCA model is generated using only the Sample data. This removes any possible bias from
Project in Blank samples	the QC, Blank, or Reference samples. The QC, Blank, or Reference sample data can
Project in Reference samples	applied to the PCA model (projected through) and plotted with the Sample data.